

aycan universal thin client viewer

Connect and Read from Anywhere

An FDA cleared medical device that can be used for primary diagnosis or review purposes, the aycan universal thin client viewer ensures a prompt and reliable way to search, view, and analyze medical images, signals, and video files on computers, smart phones, and tablets.

System Overview





HTML 5 Zero-Footprint DICOM Viewer

Runs on any HTML5 browser and doesn't require installation on the client device. The aycan application server handles the communication with the hospital systems (HIS/RIS/PACS/EMR) and does image preparation for streaming to the aycan viewer.

Supports Lightweight Directory Access Protocol (LDAP), HTTPS Security integration as well as Single Sign-On (SSO) through security tokens to allow fast and flexible usage of the system for physicians. Proprietary ways of integration to user management and access control systems is available on a project specific basis.

Key Benefits

- Connect and read from anywhere
- Supports all types of studies
- Fast responsive interface
- Native DICOM image quality
- Fast download
- Fully enabled on mobile devices
- FDA 510(k) cleared

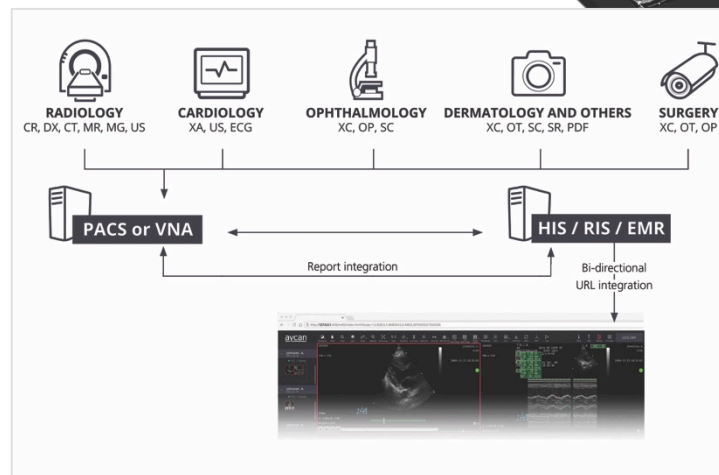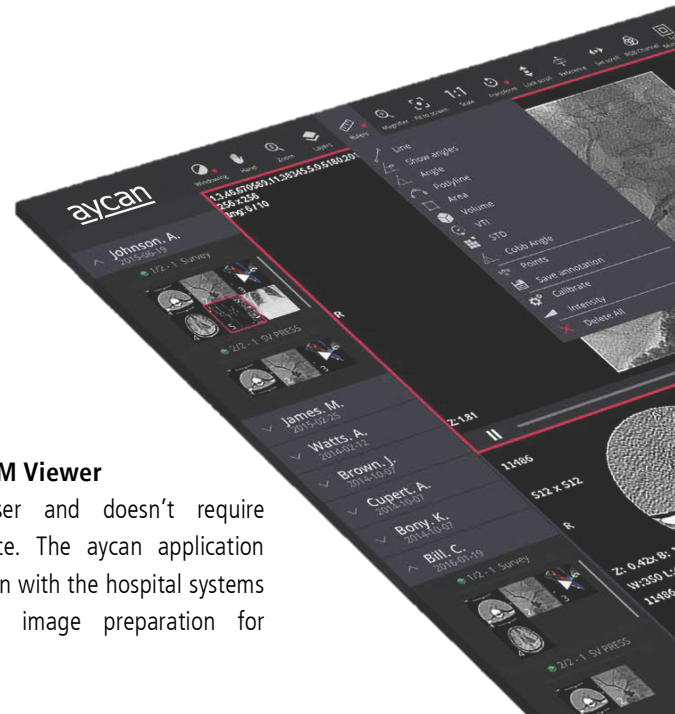

See a Live Demo/Get a Free Trial

Contact us today.
585.473.1350